

EASY GO VAC ASPIRATOR SERIES



Choice of Power Source

The EasyGoVac models give you the option of a Lithium battery for faster re-charge times and no memory, a lead acid battery or an AC unit for lightweight non-portable use.

High Flow Rates

Even with filters, all models have high flow rates of between 20 and 26 LPM, with vacuum levels of 51 – 533mmHg.

Canister Options

All units are available with a dedicated 1200cc re-useable canister or a standard 800cc canister.

Built-in Battery Charger



All battery models have a built-in charger, eliminating the need to carry a heavy external charger which can also be lost. Everything is inside the unit, to charge it, all you need to do is plug it in.

Easy to Use

All models have an easy to read “power source indicator”, “battery test indicator” (Li model only) “low battery warning light” and a “fully charged indicator” to make operation by the caregiver and/or patient easier. A convenient carry bag is included with all battery models.

***Proudly
Made in***



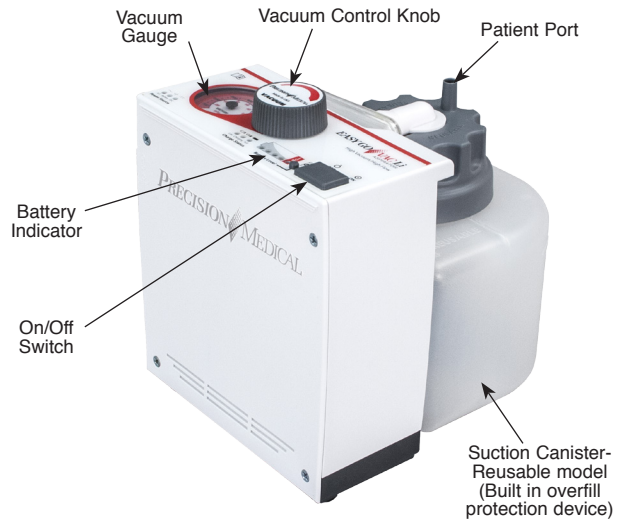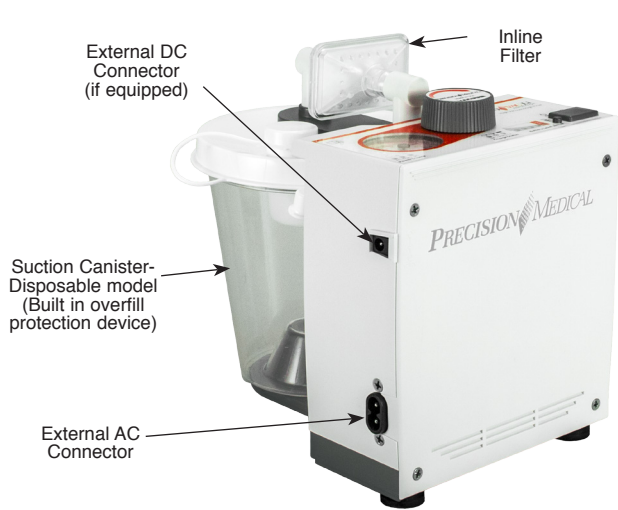

SPECIFICATIONS

		PM66AC, PM66Li		PM66S	
		Disposable Bottle (800cc)	Reusable Bottle (1200cc)	Disposable Bottle (800cc)	Reusable Bottle (1200cc)
Dimensions:	Length	8.60 in. (21.8 cm)	7.45 in. (18.9 cm)	9.43 in. (23.9 cm)	8.28 in. (21.0 cm)
	Width	6.30 in. (16.0 cm)	6.30 in. (16.0 cm)	6.30 in. (16.0 cm)	6.30 in. (16.0 cm)
	Height	8.02 in. (20.4 cm)	7.38 in. (18.7 cm)	8.02 in. (20.4 cm)	7.38 in. (18.7 cm)
Weight :	(Device)			4.90 lbs (2.22 kg)	5.00 lbs (2.27 kg)
	PM66AC:	3.25 lbs (1.47 kg)	3.35 lbs (1.52 kg)		
	PM66Li:	3.70 lbs (1.67 kg)	3.80 lbs (1.72 kg)		
	(Shipping)			7.20 lbs (3.27 kg)	7.35 lbs (3.33 kg)
	PM66AC:	4.50 lbs (2.04 kg)	4.65 lbs (2.11 kg)		
	PM66Li:	6.00 lbs (2.72 kg)	6.15 lbs (2.79 kg)		

Vacuum:
2 - 21 inHg (51 - 533 mmHg)

Flow:
≥20 LPM (Unrestricted Free Flow)

Battery Duration:
2 hours

Battery Re-Charge Time:
2 hours

AC Power Source:
100-240 VAC 50/60 Hz

Warranty:
2 years on unit
6 months on battery

Specifications are subject to change without prior notice.

ORDERING INFORMATION

Part Number	Description
PM66AC	EasyGoVac Aspirator (AC only); Carry Bag sold separately
PM66Li	EasyGoVac Aspirator with Lithium Battery includes Carry Bag
PM66S	EasyGoVac Aspirator with Sealed Lead Acid Battery includes Carry Bag

When ordering a PM66 Aspirator add:
-D for Disposable Canister
-R for Reuseable Canister

Part Number	Description
1955	Patient Suction Tubing
502519-10	800cc Canister, (Case-10 ea.)
507213	Auto DC Power Cord
506961	AC Power Cord
507263	Carry Bag
507754	Inline Filter (Bacteria)
506851	Inline Filter (Hydrophobic)
507521	800cc Inlet Adapter Kit
507522	1200cc Canister



Precision Medical, Inc.
300 Held Drive • Northampton, PA 18067 USA
Phone: (+001) 610-262-6090 Toll-Free Phone: 1-800-272-7285
Fax: (+001) 610-262-6080 Toll-Free Fax: 1-800-352-1240
www.precisionmedical.com