

Oxygen Supply Manifold



Emergency O₂ Supply When You Need it!

It's Complete...

Six oxygen masks, one 20 ft hose and the manifold supplied in a hard plastic case. Ready to use in any emergency situation requiring multiple O₂ sources.

It's Easy to Use...

Each large control knob has easy-to-read numbers and a color-coded graduated scale. Each end has a DISS male with a check valve making the connection to the oxygen supply easier.


It's Durable...

The manifold and control valves are made of durable aluminum. Rubber feet on the bottom protect the unit from the elements as well as from sliding.

It's Self Contained...

The entire unit is self-contained in the hard plastic case allowing for easy storage in the ambulance or emergency vehicle. It allows quick movement to the set-up location.



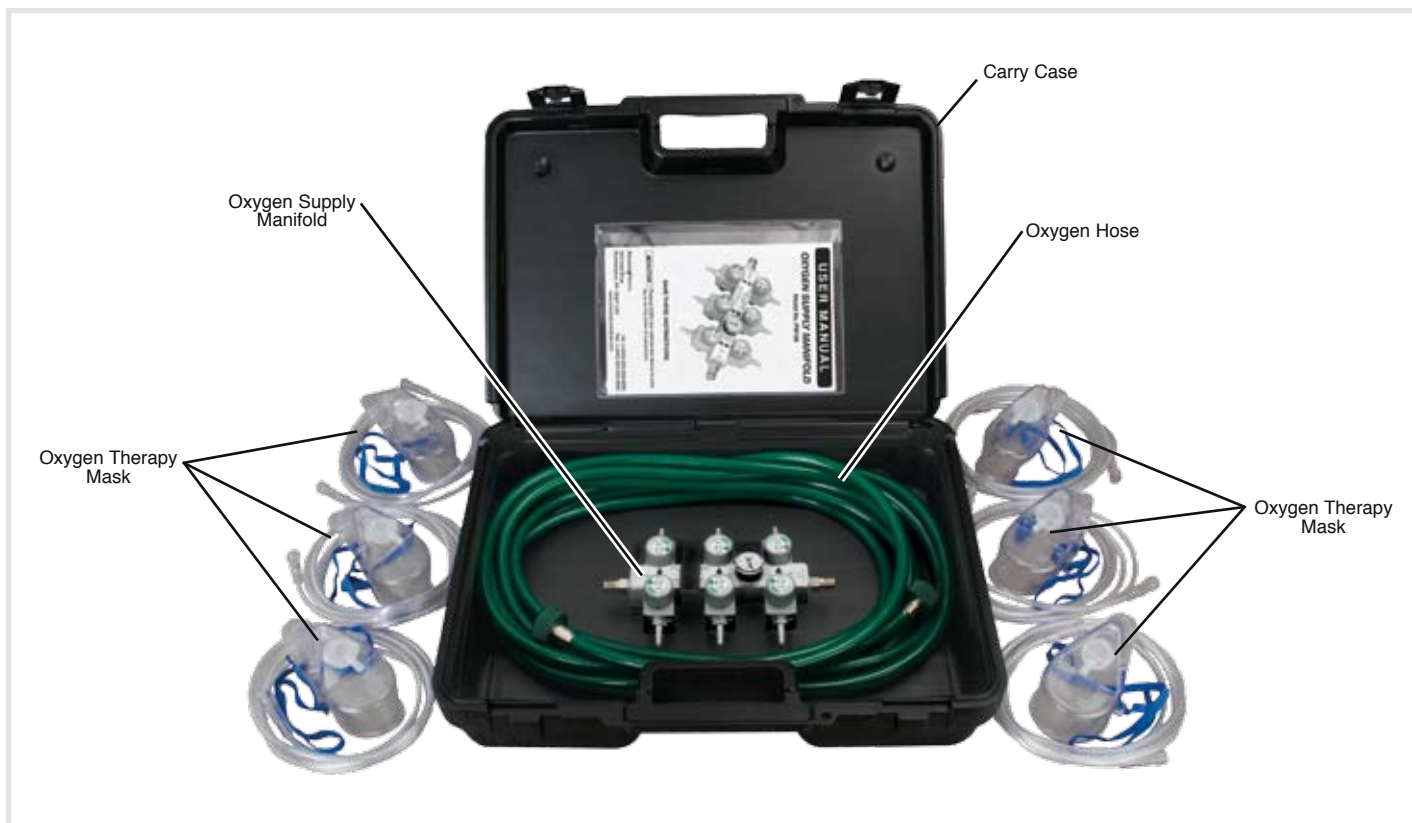

SPECIFICATIONS

Dimensions:

Oxygen Supply Manifold ONLY

Length: 7.75 in (19.69 cm)

Width: 5.00 in (12.70 cm)

Height: 2.25 in (5.72 cm)

Carry Case

Length: 12.50 in (31.75 cm)

Width: 17.00 in (43.18 cm)

Height: 4.50 in (11.43 cm)

Weight: (also shipping weight)

6.5 lbs (2.9 kg)

includes Oxygen Supply Manifold, one 20 ft Oxygen Hose, 6 Oxygen Therapy Masks and Carrying Case

Supply Pressure Requirement:

50 psi

Flow Settings:

0-15 lpm variable

ORDERING INFORMATION

Part#	Description
PM106	Oxygen Supply Manifold, one 20 ft Oxygen Hose, 6 Oxygen Therapy Masks and Carrying Case
505137-50	Oxygen Therapy Mask (case of 50)
505136	20 ft Oxygen Supply Hose, DISS Female x DISS Female



ISO 13485
Certified

Precision Medical, Inc.
300 Held Drive • Northampton, PA 18067
Phone: 610-262-6090 Toll-Free Phone: 1-800-272-7285
Fax: 610-262-6080 Toll-Free Fax: 1-800-352-1240
www.precisionmedical.com